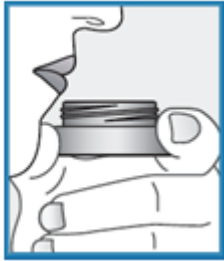

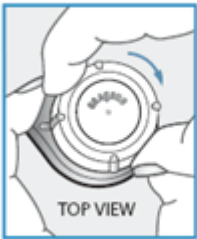

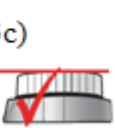



Oragene®•DNA Self-Collection Kit User Instructions (OG-250 Disc Format)

<p>1)</p> 	<p>2)</p> 	<p>3a)</p>  <p>3b)</p>  <p>3c)</p> 	<p>4)</p> 
<p>Do not eat, drink, smoke or chew gum for 30 minutes before giving your saliva sample.</p>			
<p>1 Spit saliva into the empty container (Picture #1)</p>	<p>2 Spit until the amount of liquid saliva (not bubbles) reaches the level shown in Picture #2</p>	<p>3 Put the container on a flat surface. Screw the cap onto the container. Make sure that the cap is closed tightly, as shown in picture #3c.</p> <p>If you have trouble closing the cap all the way, turn the cap slightly counter-clockwise then try to close it again.</p>	<p>4 Mix gently for at least 10 seconds.</p> <p>Recycle outer packaging.</p>

User Tips:

- Do NOT remove plastic film from the lid.
- Be patient: On average, it takes approximately 2 to 5 minutes to provide a saliva sample.
- Some people may find it hard to spit the recommended amount of saliva. To make more saliva, close your mouth and wiggle your tongue or rub your cheeks. It is easier to spit more if you place a few grains of plain white sugar on your tongue.
- Finish spitting within 30 minutes and immediately close the container.

Intended Use: This product is designed for the safe collection of DNA samples from human saliva.

Contents: The white lid contains 2 mL of Oragene•DNA liquid. Before use, the solution in the lid should be clear and colorless.

Warnings: Do not ingest the Oragene•DNA liquid. Wash with water if the Oragene•DNA liquid comes in contact with eyes or skin.

Storage: Store at room temperature 15-30°C (59-86°F).